

JOY TEST

MONO Rapid Test Cassette (Whole Blood/Serum/Plasma)

Package Insert

REF IMO-402

English

A rapid test for the diagnosis of Infectious Mononucleosis (IM) to detect Infectious Mononucleosis heterophile antibodies qualitatively in human whole blood, serum or plasma.
For professional *in vitro* diagnostic use only.

INTENDED USE

The MONO Rapid Test Cassette (Whole Blood/Serum/Plasma) is a rapid chromatographic immunoassay for the qualitative detection of Infectious Mononucleosis heterophile antibodies in human whole blood, serum or plasma as an aid in the diagnosis of Infectious Mononucleosis.

SUMMARY

Infectious Mononucleosis (IM) is caused by the Epstein-Barr virus, which is a member of the herpesvirus family. Symptoms of IM are fever, sore throat and swollen lymph glands. In very rare cases, heart or central nervous system problems may occur. Diagnosis of IM is made based on the presence of heterophile antibodies. Infectious Mononucleosis heterophile antibodies belong to the IgM class. They are present in 80-90% of acute IM cases and can be detected in 60-70% of patients during the first week of clinical illness.^{1,2,3,4}

The MONO Rapid Test Cassette (Whole Blood/Serum/Plasma) is a simple test that utilizes an extract of bovine erythrocytes to qualitatively and selectively detect Infectious Mononucleosis heterophile antibodies in whole blood, serum or plasma in minutes.

PRINCIPLE

The MONO Rapid Test Cassette (Whole Blood/Serum/Plasma) is a qualitative, lateral flow immunoassay for the detection of IM heterophile antibodies in whole blood, serum or plasma. In this test, bovine erythrocyte extracted antigen is immobilized in the test line region of the test. During testing, the specimen reacts with bovine erythrocyte extracted antigen coated particles that have been applied to the label pad. This mixture migrates chromatographically along the length of the test and interacts with the immobilized bovine erythrocyte extracted antigen. If the specimen contains IM heterophile antibodies, a colored line will appear in the test line region, indicating a positive result. If the specimen does not contain IM heterophile antibodies, a colored line will not appear in this region indicating a negative result. To serve as a procedural control, a colored line will always appear in the control line region indicating that proper volume of specimen has been added and membrane wicking has occurred.

REAGENTS

The test contains bovine erythrocyte extracted antigen-coated particles and bovine erythrocyte extracted antigen-coated membrane.

PRECAUTIONS

- For professional *in vitro* diagnostic use only. Do not use after expiration date.
- The test must remain in the sealed pouch until use.
- Do not eat, drink or smoke in the area where the specimens or kits are handled.
- Do not use test if pouch is damaged.
- Handle all specimens and controls as if they contain infectious agents. Observe established precautions against microbiological hazards throughout testing and follow standard procedures for proper disposal of specimens and controls.
- Human plasma used in the Positive and Negative Controls was tested by ELISA for the presence of antibodies to human immunodeficiency virus type HIV-1/HIV-2, as well as Hepatitis B surface antigen (HBsAg) and anti-HCV, and found to be negative. Nevertheless, caution should be used in handling and disposing of these items.
- Wear protective clothing such as laboratory coats, disposable gloves and eye protection when specimens are being tested.
- The used test should be discarded according to local regulations.
- Humidity and temperature can adversely affect results.

STORAGE AND STABILITY

Store as packaged in the sealed pouch either at room temperature or refrigerated (2-30°C). The test is stable through the expiration date printed on the sealed pouch. The test must remain in the sealed pouch until use. **DO NOT FREEZE.** Do not use beyond the expiration date.

SPECIMEN COLLECTION AND PREPARATION

- The MONO Rapid Test Cassette (Whole Blood/Serum/Plasma) can be performed using whole blood (from venipuncture or fingerstick), serum or plasma.
- To collect Venipuncture Whole Blood specimens: Collect anti-coagulated blood specimen (sodium or lithium heparin, potassium or sodium EDTA, sodium oxalate, sodium citrate) following standard laboratory procedures.
- To collect Fingerstick Whole Blood specimens:
 - Wash the patient's hand with soap and warm water or clean with an alcohol swab. Allow to dry.
 - Massage the hand without touching the puncture site by rubbing down the hand towards the fingertip of the middle or ring finger.
 - Puncture the skin with a sterile lancet. Wipe away the first sign of blood.
 - Gently rub the hand from wrist to palm to finger to form a rounded drop of blood over the puncture site.
 - Add the Fingerstick Whole Blood specimen to the test by using a capillary tube:
 - Touch the end of the capillary tube to the blood until filled to approximately 50 µL. Avoid air bubbles.
 - Place the bulb onto the top end of the capillary tube, then squeeze the bulb to dispense the whole blood to the specimen well (S) of the Test Cassette.
- Separate serum or plasma from blood as soon as possible to avoid hemolysis. Use only clear, non-hemolyzed specimens.
- Testing should be performed immediately after specimen collection. Do not leave the specimens at room temperature for prolonged periods. Serum and plasma specimens may be stored at 2-8°C for up to 3 days. For long-term storage, specimens should be kept below -20°C. Whole blood collected by venipuncture should be stored at 2-8°C if the test is to be run within 2 days of collection. Do not freeze whole blood specimens. Whole blood collected by fingerstick should be tested immediately.
- Bring specimens to room temperature prior to testing. Frozen specimens must be completely thawed and mixed well prior to testing. Specimens should not be frozen and thawed repeatedly.
- If specimens are to be shipped, they should be packed in compliance with local regulations covering the transportation of etiologic agents.

MATERIALS

- Materials Provided
- Test cassette
 - Sterile Lancet
 - Buffer
 - Alcohol Pad
 - Specimen collection containers (for venipuncture whole blood)
 - Dropper
 - Package insert

- Timer
- Centrifuge
- Heparinized capillary tubes and dispensing bulb

DIRECTIONS FOR USE

Allow the test, specimen, buffer and/or controls to equilibrate to room temperature (15-30°C) prior to testing.

- Remove the test cassette from the foil pouch and use it as soon as possible. Best results will be obtained if the assay is performed within one hour.
- Place the Test Cassette on a clean and level surface.

For Serum or Plasma specimens:

Hold the dropper vertically and transfer **1 drop of serum or plasma (approximately 25 µL)** to the specimen well (S) of the Test Cassette, and add **1 drop of buffer (approximately 55 µL)**, then start the timer. See illustration below.

For Venipuncture Whole Blood specimens:

Hold the dropper vertically and transfer **2 drops of whole blood (approximately 50 µL)** to the specimen well (S) of the Test Cassette, and add **1 drop of buffer (approximately 55 µL)**, then start the timer. See illustration below.

For Fingerstick Whole Blood specimens:

To use a capillary tube: Fill the capillary tube and transfer approximately **50 µL of fingerstick whole blood** specimen to the specimen well (S) of the Test Cassette, then add **1 drop of buffer (approximately 55 µL)** and start the timer. See illustration below.

- Wait for the colored line(s) to appear. Read results at **5 minutes**. Do not interpret the result after 10 minutes.

INTERPRETATION OF RESULTS


(Please refer to the illustration)

POSITIVE: * **Two colored lines appear.** One line should be in the control line region (C) and another line should be in the test line region (T).

***NOTE:** The intensity of the color in the test line region (T) will vary depending on the concentration of IM heterophile antibodies present in the specimen. Therefore, any shade of color in the test line region (T) should be considered positive.

NEGATIVE: **One colored line appears in the control line region (C).** No colored line appears in the test line region (T).

INVALID: **Control line fails to appear.** Insufficient specimen volume or incorrect procedural techniques are the most likely reasons for control line failure. Review the procedure and repeat the test with a new test. If the problem persists, discontinue using the test kit immediately and contact your local distributor.

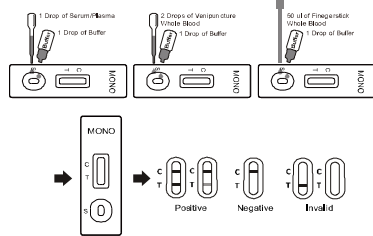

QUALITY CONTROL

A procedural control is included in the test. A colored line appearing in the control line region (C) is considered an internal procedural control. It confirms sufficient specimen volume, adequate membrane wicking and correct procedural technique.

In addition to your laboratory's standard quality control procedures, it is recommended that a positive and negative external control be tested at least once within each test kit and by each operator performing testing within a kit. This will verify that the reagents and test are working properly and the operator is able to correctly perform the test procedure. External positive and negative controls are supplied in the kit.

Procedure for External Quality Control Testing

- Holding the bottle vertically, add **1 full drop (approximately 40 µL)** of positive or negative control solution to the specimen well (S) of the Test Cassette, and add **1 drop of buffer (approximately 55 µL)**.
- Continue with Step 3 of Directions For Use.
- If the controls do not yield the expected results, do not use the test results. Repeat the test or contact your distributor.

LIMITATIONS

- The MONO Rapid Test Cassette (Whole Blood/Serum/Plasma) is for *in vitro* diagnostic use only. The test should be used for the detection of Infectious Mononucleosis antibodies in whole blood, serum or plasma specimens only. Neither the quantitative value nor the rate of increase in Infectious Mononucleosis antibody concentration can be determined by this qualitative test.
- The MONO Rapid Test Cassette (Whole Blood/Serum/Plasma) will only indicate the presence of Infectious Mononucleosis antibodies in the specimen and should not be used as the sole criteria for the diagnosis of Infectious Mononucleosis infection.
- As with all diagnostic tests, all results must be interpreted together with other clinical information available to the physician.
- If the test result is negative and clinical symptoms persist, additional testing using other clinical methods is recommended. A negative result does not at any time preclude the possibility of Infectious Mononucleosis infection.

PERFORMANCE CHARACTERISTICS

Sensitivity

The MONO Rapid Test Cassette (Whole Blood/Serum/Plasma) has been evaluated with specimens confirmed positive or negative by a leading commercial slide agglutination test. The slide agglutination test served as the reference method for the MONO Rapid Test Cassette (Whole Blood/Serum/Plasma). The result shows that the sensitivity of the MONO Rapid Test Cassette (Whole Blood/Serum/Plasma) is 97.6% relative to the slide agglutination test.

Specificity

The MONO Rapid Test Cassette (Whole Blood/Serum/Plasma) uses an antigen that is highly specific for IM antibodies in whole blood, serum or plasma. The results show that the specificity of the MONO Rapid Test Cassette (Whole Blood/Serum/Plasma) is 97.8% relative to the slide agglutination test.

MONO Rapid Test Cassette vs. Slide Agglutination

| Method | Slide Agglutination | | Total Results |
|--------------------------|---------------------|----------|---------------|
| | Positive | Negative | |
| MONO Rapid Test Cassette | Positive | 4 | 126 |
| | Negative | 176 | 179 |
| | 3 | 180 | 305 |

Relative Sensitivity: 97.6% (95%CI*: 93.1%-99.5%)*

Relative Specificity: 97.8% (95%CI*: 94.4%-99.4%)*

Overall Accuracy: 97.7% (95%CI*: 95.3%-99.1%)*

*Confidence Intervals

Precision

Intra-Assay

Within-run precision has been determined by using 3 replicates of three specimens: a negative, a low positive and a middle positive. The negative, low positive and middle positive values were correctly identified >99% of the time.

Inter-Assay

Between-run precision has been determined by 10 independent assays on the same three specimens: a negative, a low positive and a middle positive. Three different lots of the MONO Rapid Test Cassette (Whole Blood/Serum/Plasma) have been tested using negative, low positive and middle positive specimens. The specimens were correctly identified >99% of the time.

Cross-Reactivity

RF, HBsAg, HBeAg, HbCAb, HBeAb, HCV, TB, HIV and Syphilis positive specimens were tested with the MONO Rapid Test Cassette (Whole Blood/Serum/Plasma). No cross-reactivity was observed, indicating that the MONO Rapid Test Cassette (Whole Blood/Serum/Plasma) has a high degree of specificity for human antibodies to IM.

BIBLIOGRAPHY

- Hickey SM, Strasburger VC. What Every Pediatrician Should Know About Infectious Mononucleosis In Adolescents. *Pediatr Clin North Am.* 1997; 44(6):1541-56
- Omori M. Mononucleosis. 2002. <http://www.emedicine.com/EMERG/topic309.htm>
- Linde A. Diagnosis of Epstein-Barr virus-related diseases. *Scand J Infect Dis Suppl.* 1996; 100:83-8
- Papesch M, Watkins R. Epstein-Barr virus infectious mononucleosis. *Clin Otolaryngol.* 2001; 26(1): 3-8
- CDC National Center for Infectious Diseases. EBV & IM: <http://www.cdc.gov/ncidod/diseases/ebv.htm>

| | |
|--|---|
| | Caution |
| | For <i>in vitro</i> diagnostic use only |
| | Store between 2-30°C |
| | Do not use if package is damaged |
| | Authorized representative in EU |
| | Catalog # |
| | Tests per kit |
| | Use by |
| | Lot number |
| | Manufacturer |
| | Do not reuse |
| | Consult instructions for use |
| | Importer |
| | Distributor |

Hangzhou Alltest Biotech Co., Ltd.
#550, Yinhai Street
Hangzhou Economic & Technological Development Area
Hangzhou, 310018 P.R. China
www.alltests.com.cn Email: info@alltests.com.cn

MedNet EC-REP GmbH
Borkstrasse 10, 48163 Muenster, Germany
Czech Original Products s.r.o.
Koulouva 6, Praha 6, 160 00 - CZ
JOYMED.cz - ID: 08595771
Czech Original Products s.r.o.
Koulouva 6, Praha 6, 160 00 - CZ
JOYMED.cz - ID: 08595771

Number:
Revision Date: 2025-08-27